

Reset Linak CB20 & CO71 Version 1

Pedi Spa

MLW

MO1

SPX

SL


SLS

SPL

SLR

Step 1: Reset the bed

Press button A + A1 simultaneously for 10 seconds.
A beep sounds. When the beep tone ends, the reset is complete.

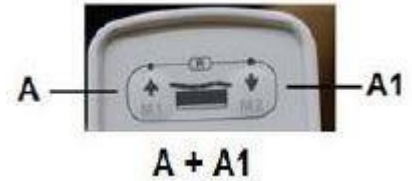


Libra Edge & Wet Table

MLX Quartz, I³Dome, Limber

Step 1: Reset the bed

Press and hold both arrow keys simultaneously for 10 seconds. A beep sounds. When the beep tone ends, the reset is complete.

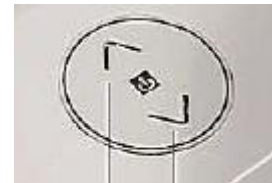

Step 2: Initialize couch

Lower both motors one after the other until couch is in end position.

MLR Handcontrol Version 1

Step 1: Reset the bed

Press and hold C1 + D1 simultaneously for 10 seconds
A beep sounds. When the beep is over, the reset is complete.

Step 2: Initialize bed

Lower C + C1 one after the other until couch is in end position.



MLR Handcontrol Version 2

Step 1: Reset the bed

Press and hold C1 + D1 simultaneously for 10 seconds
A beep sounds. When the beep is over, the reset is complete.




Step 2: Initialize bed

Lower C + C1 one after the other until couch is in end position.

PLS (Podo)


Step 1: Reset the bed

Keep buttons L1 + R2 pressed simultaneously for 10 seconds.
A beep sounds. When the beep tone ends, the reset is complete.

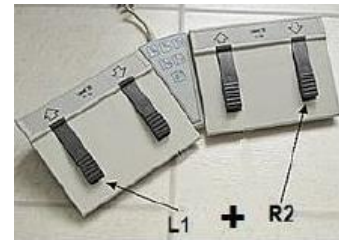

Lina Select (Alu, Wood, Static, Podo 5)

Step 1: Reset the bed

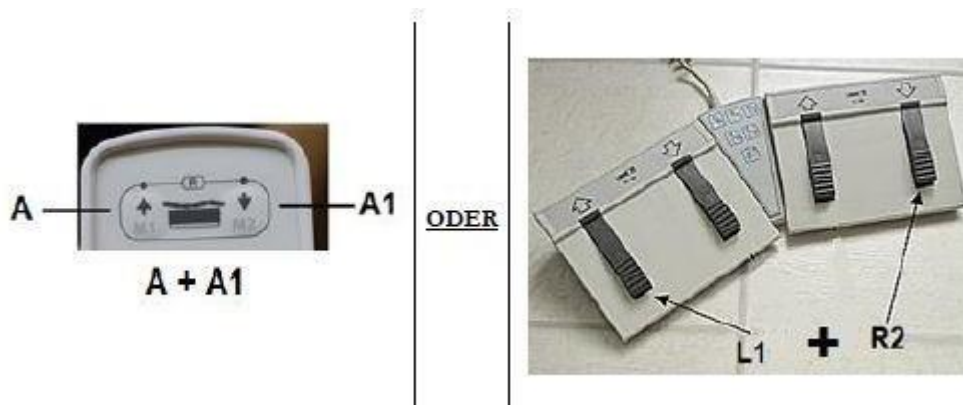

For manual operation:

Press button A + A1 simultaneously for 10 seconds.

A signal tone sounds. When the signal tone ends, the reset is complete.

By foot switch:

Press and hold buttons L1 + R2 simultaneously for 10 seconds.

A signal tone sounds. When the signal tone ends, the reset is complete.


After resetting the couch, initialization is necessary. To do this, please extend all motors completely once and retract them again.

For the initialization of the lifting column, please follow the instructions for the corresponding couches.

Reset Linak CO71 Version 2 with Double-Click System and built-in heating control

Step 1:


Press buttons T 11 + T 12 on the hand control simultaneously and keep them pressed until the LED (diode) flashes three times.

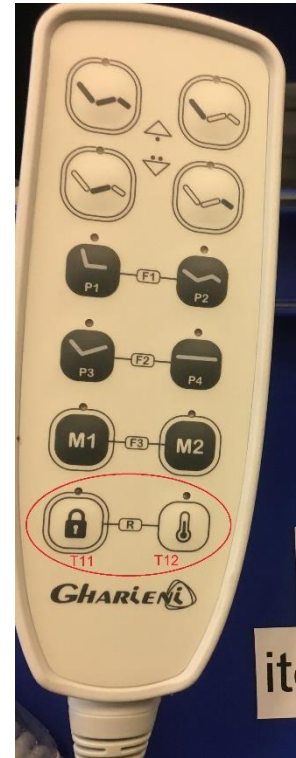

Reset TiMotion

Step 1:

Press buttons T 11 + T 12 on the hand control simultaneously and keep them pressed until the LED (diode) flashes three times.

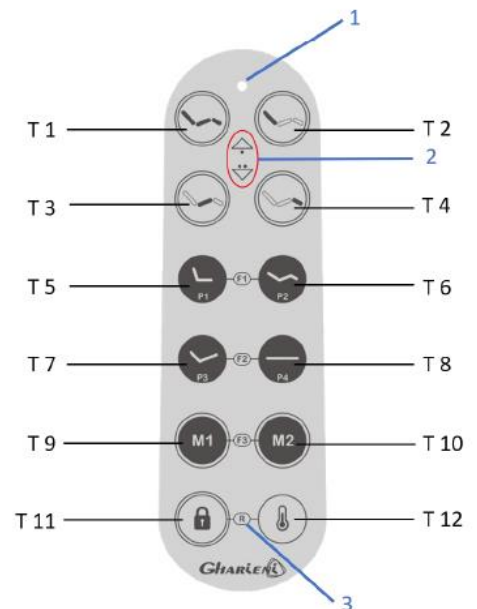

After resetting the couch, initialization is necessary. To do this, please extend all motors completely once and retract them again.

For the initialization of the lifting column, please follow the instructions for the corresponding couches.